



DRSS Received
June 20, 2025

ADAMS Accession # ML25175A081

June 16, 2025

U.S. Nuclear Regulatory Commission
Attn: Document Control Desk
Washington, DC 20555-001

RE: Reply to Notice of Violation, NRC Inspection Report 030-39297/2024-001

To Whom It May Concern:

Please find below Advanced Isotopes of Montana's response to the Notice of Violation issued on May 20, 2025.

The Notice of Violation identified two (2) Severity Level IV violations: Inadequate policies and procedures regarding Ge-68/Ga-68 generator usage to prepare Ga-68 radiopharmaceuticals and failure to ensure that hazmat employees are tested on training regarding DOT requirements. Detailed descriptions and remedies can be found in supplemental documents.

Advanced Isotopes of Montana is serious about addressing and correcting the issues identified by the NRC inspection report and improving our program. We appreciate the thoroughness and professionalism demonstrated by the NRC inspectors during their visit to our pharmacy.

Please contact me at 208-339-5120 if you have any questions.

Sincerely,

A handwritten signature in black ink, appearing to read 'Michael Hart', is written over a light gray rectangular background.

Michael Hart
Radiation Safety Officer
Advanced Isotopes of Montana, LLC



Supplements:

Supplement 1 – Reply to Violation: Ge-68/Ga-68 Generator Procedures

Supplement 2 – Reply to Violation: DOT/HAZMAT Training Testing

Attachments:

Attachment 1 – Updated Ge-68/Ga-68 Procedures

Attachment 2 – Updated Training Certificates



Supplement 1
Reply to Violation: Ge-68/Ga-68 Generator Procedures

NRC Statement of Violation (Pulled directly from the Notice of Violation)

Contrary to the above, from October 24, 2023, to December 10, 2024, the licensee failed to conduct its program in accordance with the statements, representations, and procedures contained in the documents, including any enclosures, listed in the license. Specifically, the licensee's "Policies and Procedures," Section 2.6, "Ge-68/Ga-68 Generator Elution," failed to address licensee commitments for the use of Eckert & Ziegler (EZ) GalliaPharm Ge-68/Ga-68 generators to prepare Ga-68 radiopharmaceuticals, and failed to address the EZ GalliaPharm "Instructions for Use." The licensee's policies and procedures failed to contain adequate guidance or instruction on:

- (1) Following the EZ GalliaPharm "Instructions for Use," including: unpacking the generator, optimal generator positioning, generator preparation including attachment of specified inlet and outlet port lines with proper manifolds and eluent, first elution (conditioning procedure), continuous routine elution, and exchange of sterile hydrochloric acid container.*
- (2) Eluting the generator in accordance with the manufacturer's stated frequency, including decay correcting the measured activity to the starting time of the elution and taking steps to minimize the concentration of Ge-68 in the eluate.*
- (3) Not using an expired generator.*
- (4) Following the manufacturer's instructions for "Shelf Life of a GalliaPharm Generator."*
- (5) Breakthrough testing and frequency, including appropriate methods for the determination of Ge-68 breakthrough that will detect whether the eluate exceeds the EZ GalliaPharm generator recommended breakthrough limit of not more than 0.001 percent of the eluted Ga-68 activity.*
- (6) Maintaining records of the breakthrough tests that include the ratio of the measured activity of Ge-68 per Ga-68, corrected for the time of elution, the time and date of the elution, the time and date of the measurement, and the name of the individual who made the measurement.*
- (7) Notifying by telephone the NRC Operations Center and EZ within 7 calendar days after discovery of a GalliaPharm generator that is unable to meet the manufacturer's stated Ge-68 breakthrough limits.*
- (8) Sending a written report to the appropriate NRC Regional Office within 30 days after discovery of a generator that is unable to meet the manufacturer's stated breakthrough limits of Ge-68.*
- (9) Emergency procedures for leaking or damaged generators.*

**Reason for the Violation**

Advanced Isotopes of Montana agrees that this was a violation of NRC requirements and commitments made by Advanced Isotopes of Montana.

The reason for the violation was due to the storage of the policies and procedures in an email system that had its services revoked and the use of outdated policies and procedures that did not take into account the Ge-68/Ga-68 generators.

Corrective steps that have been taken and results achieved

Policies and procedures were updated to include adequate guidance and instruction on the 9 items identified in the notice of violation. The applicable policies and procedures are included in this submission.

Corrective actions that will be taken to avoid further violations

The policies and procedures will be overhauled and reviewed for accurateness. Policies and procedures will be reviewed by management and full-time staff annually, or when new products are added to Advanced Isotopes of Montana's formulary, to ensure that practices reflect procedures and to identify areas of improvement.

Date when corrective action will be completed

Updated policies and procedures went into effect on 6/16/2025.



Supplement 2
Reply to Violation: DOT/HAZMAT Training Testing

NRC Statement of Violation (Pulled directly from the Notice of Violation)

10 CFR 71.5(a) requires, in part, that each licensee who transports licensed material on public highways shall comply with the applicable requirements of the U.S. Department of Transportation (DOT) regulations in 49 CFR Parts 171 through 180 appropriate to the mode of transport.

49 CFR 172.702(d) requires that a hazmat employer shall ensure that each of its hazmat employees is tested by appropriate means on the training subjects covered in 49 CFR 172.704.

Contrary to the above, from November 12, 2024, to December 10, 2024, the licensee transported licensed material on public highways and failed to ensure that each of its hazmat employees is tested by appropriate means on the training subjects covered in 49 CFR 172.704. Specifically, the licensee failed to ensure that a hazmat employee was tested on the training subjects covered in 49 CFR 172.704. The hazmat employee independently performed job functions involving the transportation of Class 7 (radioactive) materials on public highways in areas of NRC jurisdiction.

Reason for the Violation

Advanced Isotopes of Montana agrees that this was a violation of DOT/NRC regulations. The individual identified in the citation had received training, but we had lost the exam proving completion of the course.

Corrective steps that have been taken and results achieved

The individual identified in the citation completed a new DOT training test for 2024. In 2025, the individual completed DOT/Hazmat training provided through Nuclear Education Online and Curium Pharma LLC.

Corrective actions that will be taken to avoid further violations

Future exams will be completed and scanned into the employee's file prior to the employee being allowed to deliver without appropriate supervision.

Date when corrective action will be completed

Corrective action was completed in February 2025 by all employees.

2.6. Ge-68/Ga-68 Generator Elution

2.6.1. Purpose

- 2.6.1.1. This procedure defines the steps to elute a Germanium-68/Gallium-68 generator. The generator eluate is the starting material for radiopharmaceutical compounding of the Gallium-68 based imaging agents.

2.6.2. Rationale

- 2.6.2.1. Proper generator elution allows the pharmacist to work with eluates of desired activities and concentrations in preparing kits and dispensing doses.

2.6.3. Equipment

- | | | |
|----------|-----------------------|----------------------------|
| 2.6.3.1. | Ge-68/Ga-68 generator | Lead shielded elution vial |
| | Dose calibrator | Remote handling device |

2.6.4. Procedure

- 2.6.4.1. After receipt of the generator and performance of appropriate radioactive material receiving procedures (wipes, surveys, log-in, etc.), remove generator from its shipping carton and place it in a secondary shield. The secondary shield is usually provided by the generator manufacturer and should be in place before receiving the generator.
- 2.6.4.1.1. If the shipping container shows signs of damage or droppage, check for visible liquid on the generator and visible damage. Gently tilt the generator 90 degrees and listen for moving internal parts, if internal noise is heard, the generator may be damaged.
- 2.6.4.1.2. If a generator is suspected or determined to be leaking or damaged, **DO NOT USE**.
- 2.6.4.1.2.1. Call RSO immediately if the generator is damaged or leaking.
- 2.6.4.1.2.2. Follow emergency spill precautions in the Radiation Safety Manual (also posted in the Hot Lab).
- 2.6.4.1.2.3. The RSO will call the manufacturer and contain any spills.
- 2.6.4.2. Following the instructions given with the sterile HCl lines, aseptically install the input and output lines for elution.
- 2.6.4.3. Spike a new 0.1 mol/L HCl eluent bag
- 2.6.4.3.1. Each bag should allow for ~40 elutions
- 2.6.4.3.2. When changing bags, ensure there is no air making its way into the inlet side of the generator
- 2.6.4.4. Prior to elution, ensure that the generator has been eluted in the last 48-72 hours.
- 2.6.4.4.1. If it has not, the generator must be eluted using 10 mL to ensure that there is no Ge-68 breakthrough.

2.6.4.5. NOTE: The generator is good for 365 days after calibration date or 1650 elutions

2.6.4.5.1. DO NOT USE EXPIRED ELUTIONS FOR PATIENT DOSES

- 2.6.4.6. Calculate routine Gallium-68 utilization requirements. Record date of receipt of generators to be eluted.
- 2.6.4.7. Determine the desired eluate volume (5 mL for Illucix, 10 mL for breakthrough assay, 5-10 mL to flush every 48-72 hours.)
- 2.6.4.8. Place the elution vial in a lead shielded container.
- 2.6.4.9. Aseptically connect the shielded vial to the appropriate generator and elute for the desired time or volume.
 - 2.6.4.9.1. Note the time that the elution began and enter that time as the initial elution time.
- 2.6.4.10. Remove the vial.
- 2.6.4.11. In a properly calibrated dose calibrator, determine the total activity and record this activity and volume of each eluate in the computer.
 - 2.6.4.11.1. To determine elution activity, use the decay ratio chart laminated in the hot lab, divide by the appropriate ratio given the amount of time that has passed.
 - 2.6.4.11.2. Enter the elution start time and the calculated activity as the elution amount.
 - 2.6.4.11.3. Follow the appropriate procedure for either dose preparation, Ge-68 breakthrough assay, or dispose of vial as a waste elution.

2.6.5. Elution Ge-68 Breakthrough Assay Procedure

- 2.6.5.1. Elute the generator in 10 mL as described above.
- 2.6.5.2. Store the breakthrough elution in a shielding area for 48-72 hours before measurement.
- 2.6.5.3. Prior to manipulating the elution, record a background measurement on the multi-channel analyzer/well counter for 60 seconds using settings appropriate for reading Ge-68/Ga-68.
- 2.6.5.4. Remove an appropriate volume of elution for the vial size being used in the multi-channel analyzer/well counter and place it into a scintillation vial, centrifuge tube, or similarly sized container for measurement.
 - 2.6.5.4.1. Appropriate volumes are defined as roughly half of the volume maximum of the container being used in the multi-channel analyzer/wellcounter.
- 2.6.5.5. Measure the elution sample in the MCA/well counter using the same settings as the previous background reading.
 - 2.6.5.5.1. Counts should be at or around background levels.
- 2.6.5.6. In BioRx, go into the elutions list in processing and sort for the breakthrough elution and double click the elution.

- 2.6.5.7. Using the Convert from DPM calculation, enter the background, measurement, window, threshold, and efficiency for the reading.
 - 2.6.5.7.1. BioRx will then convert counts to micro-curies and automatically enter the value in the breakthrough calculation.
- 2.6.5.8. If the measurement is appropriate, BioRx will highlight the reading in white.
- 2.6.5.9. If the measurement exceeds 0.001%, BioRx will show an error and mark the elution as red.
 - 2.6.5.9.1. If this happens, contact the RSO ASAP.
 - 2.6.5.9.2. The RSO will notify the NRC and generator manufacturer within 7 calendar days of the breakthrough failing.
 - 2.6.5.9.2.1. A written report will be submitted within 30 days.
- 2.6.5.10. Breakthrough results will be kept in BioRx for review for at least 3 years.

9. Emergency Procedures

A copy of these procedures will be posted in each area where radioactive material is used or stored. A decontamination kit is located in each room within the restricted area. The contents of this kit include: disposable gloves, housekeeping gloves, disposable lab gowns, disposable shoe covers, absorbent paper with plastic backing, masking tape, plastic trash bags with twist ties, "Radioactive Material" labeling tape, marking pen, pre-strung "Radioactive Material" labeling tags, box of wipes, instructions for "Emergency Procedures," clipboard with a copy of the Radioactive Spill/Incident Report Form, and a pen or pencil.

9.1. Definition of Major vs Minor Spills

9.1.1. All major spills are defined as any activity equal to or above the estimated values in the table below, including leaking or damaged Mo-99 or Ge-68 generators.

9.1.2. All minor spills are defined as any activity below the estimated values in the table below

9.1.3.

Radionuclide	Millicuries (estimate)	Radionuclide	Millicuries (estimate)
Fluorine-18	10	Iodine-123	10
Gallium-67	100	Iodine-131	1
Indium-111	10	Gallium-68	10
Yttrium-90	100	Copper-64	10
Technetium-99m	100	Thallium-201	100

9.2. Minor Spills

9.2.1. Notify persons in the area that the spill has occurred.

9.2.2. Prevent the spread of contamination by covering the spill with absorbent paper. Paper should be dampened if solids are spilled.

9.2.3. Clean up the spill using disposable gloves, absorbent paper and remote handling tongs. Carefully fold the absorbent paper. Insert absorbent paper and all other contaminated materials, such as disposable gloves, into a plastic bag. Put the plastic bag into the appropriate radioactive waste container.

9.2.4. Resurvey the area. Check the area around the spill, hands, and clothing for contamination.

9.2.5. Report the incident to the RSO promptly.

9.2.6. The RSO will:

9.2.6.1. Follow up on the decontamination activities and document the results.

9.2.6.2. As appropriate, determine cause and corrective actions needed; consider bioassays if licensed material may have been ingested or inhaled.

9.2.6.3. If necessary, notify NRC.

9.3. Major Spills

- 9.3.1. Instructions to Workers: Clear the area. If appropriate, survey all persons not involved in the spill and vacate the room.
- 9.3.2. Prevent the spread of contamination by covering the spill with absorbent paper (paper should be dampened if solids are spilled), but do not attempt to clean it up. To prevent the spread of contamination, limit the movement of all personnel who may be contaminated.
- 9.3.3. Shield the source only if it can be done without further contamination or significant increase in radiation exposure.
- 9.3.4. Close the room and lock the doors to prevent entry. Post a sign to warn anyone trying to enter that a spill of radioactive material has occurred.
- 9.3.5. Call the RSO immediately.
- 9.3.6. Survey all personnel who could possibly have been contaminated. Decontaminate personnel by removing contaminated clothing and flushing contaminated skin with lukewarm water and then washing with a mild soap.
- 9.3.7. Allow no one to return to work in the area unless approved by the RSO.
- 9.3.8. Follow the instructions of the RSO (e.g., decontamination techniques, surveys, provision of bioassay samples, requested documentation).
- 9.3.9. The RSO will:
 - 9.3.9.1. Confirm decontamination of personnel. If decontamination of personnel was not fully successful, consider inducing perspiration by covering the area with plastic. Then wash the affected area again to remove any contamination that was released by the perspiration.
 - 9.3.9.2. Skin contamination must be evaluated to determine potential exposures. Beta-emitting radionuclides have a high potential for resulting in shallow-dose exposures in excess of regulatory limits from small (microcurie) quantities of contamination.
 - 9.3.9.3. Supervise decontamination activities and document the results. Documentation should include location and results of surveys and decontamination results.
 - 9.3.9.4. Determine root cause and needed corrective actions; consider need for bioassays if licensed material may have been ingested, inhaled, or absorbed.
 - 9.3.9.5. If necessary, notify NRC.

9.4. Fire

9.4.1. Minor Fires

- 9.4.1.1. Instructions to Workers:

- 9.4.1.1.1. If possible, immediately attempt to put out the fire by approved methods (e.g., fire extinguisher) if other fire hazards or radiation hazards are not present.
- 9.4.1.1.2. Notify all persons present to vacate the area and have one individual immediately call the RSO and fire department or 911 (as instructed by the RSO).
- 9.4.1.1.3. Once the fire is out, isolate the area to prevent the spread of possible contamination.
- 9.4.1.1.4. Ensure injured personnel receive medical attention. Survey all persons involved in combating the fire for possible contamination.
- 9.4.1.1.5. Decontaminate personnel by removing contaminated clothing and flushing contaminated skin with lukewarm water, then washing with a mild soap.
- 9.4.1.1.6. In consultation with the RSO, determine a plan of decontamination and the types of protective devices and survey equipment that will be necessary to decontaminate the area.
- 9.4.1.1.7. Allow no one to return to work in the area unless approved by the RSO.
- 9.4.1.1.8. Follow the instructions of the RSO (e.g., decontamination techniques, surveys, provision of bioassay samples, requested documentation).
- 9.4.1.2. The RSO will:
 - 9.4.1.2.1. Notify emergency medical personnel of any injured individuals who may be contaminated. Provide radiation safety assistance (e.g., monitoring) as needed or requested.
 - 9.4.1.2.2. Supervise decontamination activities at the facility. If decontamination of personnel was not fully successful, consider inducing perspiration by covering the area with plastic. Then wash the affected area again to remove any contamination that was released by the perspiration.
 - 9.4.1.2.3. Consult with fire safety officials to ensure that there is no likelihood of fire restarting and that it is safe to re-enter the building.
 - 9.4.1.2.4. Determine cause and needed corrective actions; consider need for bioassays if licensed material may have been ingested or inhaled. Document the incident.
 - 9.4.1.2.5. If necessary, notify NRC.
- 9.4.2. Fires, Explosions, or Major Emergencies
 - 9.4.2.1. Instructions to Workers:
 - 9.4.2.1.1. Notify all persons in the area to leave immediately.
 - 9.4.2.1.2. Notify the fire department or 911.

- 9.4.2.1.3. Notify the RSO and other facility safety personnel.
- 9.4.2.1.4. Ensure injured personnel receive medical attention.
- 9.4.2.1.5. Upon arrival of firefighters, inform them where radioactive materials are stored or where radionuclides were being used; inform them of the present location of the licensed material and the best possible entrance route to the radiation area, as well as any precautions to avoid exposure or risk of creating radioactive contamination by use of high pressure water, etc.
- 9.4.2.1.6. Allow no one to return to work in the area unless approved by the RSO.
- 9.4.2.1.7. Follow the instructions of the RSO (e.g., decontamination techniques, surveys, provision of bioassay samples, requested documentation).
- 9.4.2.2. The RSO will:
 - 9.4.2.2.1. Notify emergency medical personnel of any injured individuals who may be contaminated.
 - 9.4.2.2.2. Provide radiation safety assistance (e.g., monitoring) as needed or requested.
 - 9.4.2.2.3. Coordinate activities with local fire departments or other emergency personnel.
 - 9.4.2.2.4. Consult with the firefighting personnel or other emergency personnel and set up a controlled area where personnel can be surveyed for contamination of their protective clothing and equipment after the fire is extinguished.
 - 9.4.2.2.5. Once the fire is extinguished, provide assistance to firefighters or other emergency personnel who may need to re-enter restricted areas to determine the extent of the damage to the licensed material use and storage areas.
 - 9.4.2.2.6. To the extent practical, assist firefighters and emergency personnel in maintaining their exposures ALARA if the fire resulted in a significant release of radioactive material or loss of shielding capability, such that excessive radiation levels (greater than 100 millirems per hour) are created.
 - 9.4.2.2.7. Perform thorough contamination surveys of firefighters and emergency personnel and their equipment before they leave the controlled area, and decontaminate if necessary.
 - 9.4.2.2.8. Supervise decontamination activities.
 - 9.4.2.2.9. Consider bioassays if licensed material may have been ingested or inhaled. Document incident.
 - 9.4.2.2.10. If necessary, notify NRC.

- 9.4.3. Copies of emergency procedures should be provided to all users. A current copy of the emergency procedures should be posted in each area where radioactive material is used.
- 9.5. Theft
 - 9.5.1. A pharmacy vehicle containing radioactive materials shall not be left unattended. If the vehicle is left unattended, then the radioactive materials shall be secured and the vehicle locked.
 - 9.5.2. If a vehicle containing radioactive material is stolen, immediately contact the RSO, the law enforcement agencies, and the NRC.
- 9.6. Motor Vehicle Accident
 - 9.6.1. A copy of "Procedures to Follow for a Motor Vehicle Accident" will be kept in all vehicles at all times.
 - 9.6.2. Standard Procedure for a Motor Vehicle Accident:
 - 9.6.2.1. Using a GM meter, determine if there is a spill or if there is contamination outside the transportation container.
 - 9.6.2.2. Contact the RSO
 - 9.6.2.3. Contact the NRC
 - 9.6.3. Begin decontamination procedures as instructed by the RSO and/or the NRC.
 - 9.6.4. Follow usual procedures as outlined for notification of police and insurance companies.

Certificate of Completion

This certifies that

jayden miles

has completed the

CURIUM™

DOT Training for AIR Shippers & Transporters of Radioactive Material/Radiopharmaceuticals
(AIR/IATA) Course

January 31, 2025

Modules covered in this course:

- Regulation of the Transport of Radioactive Materials
- Packaging and Labeling
- Documentation

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

jayden miles

has completed the

CURIUM™


DOT Training for Highway Shippers & Transporters of Radioactive Material/Radiopharmaceuticals
Radiopharmacy/Manufacturing Facility Course

January 30, 2025

Modules covered in this course:

- General Awareness of Hazardous (Radioactive) Materials
- Safety Training for Hazardous (Radioactive) Materials
- Security Awareness of Radioactive Materials
- Transporting & Delivering packages containing Radioactive Materials
- Requirement for Function-Specific Training
- Packing & Wipe Testing a Delivery Case
- Labeling Packages Containing Radioactive Materials
- Preparing Shipping Papers for Radioactive Shipments for Ground Delivery

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.

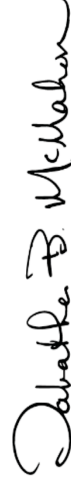

Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

Kerry clay

has completed the

CURIUM™

DOT Training for Highway Shippers & Transporters of Radioactive Material/Radiopharmaceuticals
Radiopharmacy/Manufacturing Facility Course

February 14, 2025

Modules covered in this course:

- General Awareness of Hazardous (Radioactive) Materials
- Safety Training for Hazardous (Radioactive) Materials
- Security Awareness of Radioactive Materials
- Transporting & Delivering packages containing Radioactive Materials
- Requirement for Function-Specific Training
- Packing & Wipe Testing a Delivery Case
- Labeling Packages Containing Radioactive Materials
- Preparing Shipping Papers for Radioactive Shipments for Ground Delivery

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

Kerry clay

has completed the

CURIUMTM

DOT Training for AIR Shippers & Transporters of Radioactive Material/Radiopharmaceuticals

(AIR/IATA) Course

February 20, 2025

Modules covered in this course:

- Regulation of the Transport of Radioactive Materials
- Packaging and Labeling
- Documentation

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

Michael Durnam

has completed the

CURIUM™

DOT Training for AIR Shippers & Transporters of Radioactive Material/Radiopharmaceuticals

(AIR/IATA) Course

February 20, 2025

Modules covered in this course:

- Regulation of the Transport of Radioactive Materials
- Packaging and Labeling
- Documentation

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

Michael Durnam

has completed the

CURIUM™

DOT Training for Highway Shippers & Transporters of Radioactive Material/Radiopharmaceuticals
Radiopharmacy/Manufacturing Facility Course

February 16, 2025

Modules covered in this course:

- General Awareness of Hazardous (Radioactive) Materials
- Safety Training for Hazardous (Radioactive) Materials
- Security Awareness of Radioactive Materials
- Transporting & Delivering packages containing Radioactive Materials
- Requirement for Function-Specific Training
- Packing & Wipe Testing a Delivery Case
- Labeling Packages Containing Radioactive Materials
- Preparing Shipping Papers for Radioactive Shipments for Ground Delivery

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico

Certificate of Completion

This certifies that

Michael Hart

has completed the

CURIUMTM

DOT Training for AIR Shippers & Transporters of Radioactive Material/Radiopharmaceuticals
(AIR/IATA) Course

February 14, 2025

Modules covered in this course:

- Regulation of the Transport of Radioactive Materials
- Packaging and Labeling
- Documentation

Online training has been successfully completed and tested in accordance with 49 CFR Part 172, Subpart H to prepare and transport packages of radioactive materials/radiopharmaceuticals for use in or incident to research, medical diagnosis or treatment.



Authorized Signature of the Provider
University of New Mexico